
# memmert


## ICH750L

Unparalleled homogeneity of temperature and humidity. For stability testing according to ICH - WHO - EMA - ASEAN - GMP - GLP – GCCP.



On this page, you can find all the essential technical data on the Memmert climate chamber ICH. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at <u>sales@memmert.com</u>.

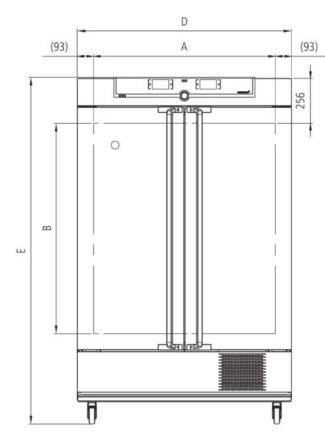

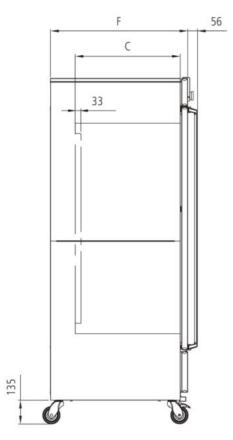

## Temperature

Working temperature range	without humidity: 0 to +60 °C
Setting temperature range	0 to +60 °C
Working temperature range	from +10°C up to +60°C with humidity and/or light
Setting accuracy temperature	0.1 °C
Temperature sensor	2 Pt100 sensors DIN Class A in 4-wire-circuit for mutual monitoring, taking over functions in case of an error

## Humidity

Humidity	humidity supply with distilled water from external tank by self-priming pump
Humidification	humidification by hot steam generator
Humidity	active humidifying and de-humidifying adjustable from 10-80 % rh with digital display of relative humidity - resolution of display 0.1 %, setting accuracy 0.5 %

## Control of standard components

Light	daylight and UV-light can be programmed individually via controller
Light	cold white light +/- 6 %, for UV-light +/- 10 %
Light	illumination unit consisting of 4 fluorescent lights with cold white light (standard illuminant D65, 6.500K) and 2 UV-lamps in the spectral range of 320-400 nm

Control technology		
ControlCOCKPIT	TwinDISPLAY. Adaptive multifunctional digital PID-microprocessor controller with 2 high-definition TFT-colour displays.	
Language setting	German, English, Spanish, French, Polish, Czech, Hungarian	
Timer	Digital backwards counter with target time setting, adjustable from 1 minute to 99 days	
Function SetpointWAIT	the process time does not start until the set temperature is reached	
Calibration	three freely selectable values each, temperature/humidity, for two freely selectable values CO2	
adjustable parameters	temperature (Celsius or Fahrenheit), fan speed, light, programme time, time zones, summertime/wintertime	

#### Ventilation

Fan

forced air circulation in working chamber adjustable from 10 to 100 % in 10 % steps

#### Communication

	programme stored in case of power failure	
Programming AtmoCONTROL software on a USB stick for programming, managing and t via Ethernet interface or USB port	d transferring programmes	

## Safety

Temperature control	mechanical temperature limiter TB, protection class 1 according to DIN 12880 to switch off the heating approx. 10°C above nominal temperature
Temperature control	over- and undertemperature monitor TWW, protection class 3.3 or adjustable temperature limiter TWB, protection class 2, selectable on display
AutoSAFETY	additionally integrated over- and undertemperature protection "ASF", automatically following the setpoint value at a preset tolerance range, alarm in case of over- or undertemperature, heating is switched off in case of overtemperature, compressor in case of undertemperature
Autodiagnostic system	integral fault diagnostics for temperature and humidity control
Alarm	visual and acoustic

## Heating concept

Air jacket	air jacket heating system for gentle all-around heating	
Heating concept	no drying-up of the load caused by the cooling device due to separation of thermal jacket from interior	
Cooling	energy-saving use of CFC-free cooling/heating system (refrigerant R134a)	
Defrosting	highly efficient automatic defrosting system	
U		

## Standard equipment

Works calibration certificate	standard values +10 °C and +37 °C, 60 % rh at +30 °C	
Door	inner glass doors	
Door	fully insulated stainless steel doors with2-point locking (compression door lock)	
Entry port	Entry port (silicone), 40 mm clear diameter, moisture tight, can be closedby a silicone stopper, standard position at the back	
Internals	2 stainless steel grid(s), electropolished	
Standard accessories	Water tank including connection hose	

## **Stainless steel interior**

Dimensions	w <sub>(A)</sub> x h <sub>(B)</sub> x d <sub>(C)</sub> : 1040 x 1200 x 600 mm (d less 33 mm for fan)
Volume	749
Max. number of internals	14
Max. loading of chamber	200 kg
Max. loading per internal	30 kg

Dimensions	w <sub>(D)</sub> x h <sub>(E)</sub> x d <sub>(F)</sub> : 1224 x 1950 x 785 mm (d +56mm door handle)
Installation	on lockable castors
Housing	rear zinc-plated steel

## Textured stainless steel casing

#### **Electrical data**

Voltage Electrical load	230 V, 50/60 Hz approx. 1550 W	
Voltage Electrical load	115 V, 50/60 Hz approx. 1550 W	

## **Ambient conditions**


Set Up	The distance between the wall and the rear of the appliance must be at least 15 cm. The clearance from the ceiling must not be less than 20 cm and the side clearance from walls or nearby appliances must not be less than 5 cm.	
Ambient temperature	15 °C to 28 °C (up to 34 °C with limited temperaturerange)	
Humidity rh	max. 70 %, non-condensing	
Altitude of installation	max. 2,000 m above sea level	
Overvoltage category	II	
Pollution degree	2	

## Packing/shipping data

Transport information	The appliances must be transported upright
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-RegNo.	DE 66812464
Dimensions approx incl. carton	w x h x d: 1330 x 2150 x 1050 mm
Net weight	approx. 249 kg
Gross weight carton	approx. 319 kg

#### Temperature-humidity working range ICH

Not all climate chambers are the same. The humidity content of the chamber load, the ambient conditions and the respective temperature-humidity working range are decisive factors in the selection of the right appliance. In the adjacent diagram, you can see the possible temperature/humidity combinations for our climate chamber ICH.Within the respective temperature-humidity range, condensation-free permanent operation is possible. To which extent condensation may occur in the threshold range depends on the humidity content of the chamber load and the ambient conditions.

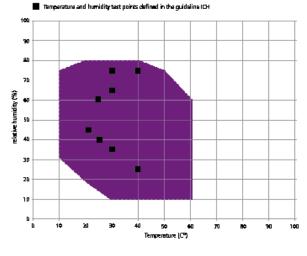

Temperature-humidity working range ICHeco/ICH

#### Standard units are safety-approved and bear the test marks

C E E H [